

#### US006042266A

# United States Patent [19]

## Cheslock et al.

# [11] **Patent Number:** 6,042,266

# [45] Date of Patent: \*Mar. 28, 2000

# [54] TYMPANIC THERMOMETER PROBE COVER

[75] Inventors: Edward P. Cheslock, Lincoln

University; Eric L. Canfield, Chester Springs; Richard K. Harris, West

Chester, all of Pa.

[73] Assignee: TruTek, Inc., West Chester, Pa.

[\*] Notice: This patent is subject to a terminal dis-

claimer.

[21] Appl. No.: 09/049,926

[22] Filed: Mar. 30, 1998

# Related U.S. Application Data

[62] Division of application No. 08/747,423, Nov. 12, 1996, Pat. No. 5,833,367.

[51] **Int. Cl.**<sup>7</sup> ...... **G01K 1/08**; A61B 6/00

[52] **U.S. Cl.** ...... **374/158**; 374/121; 600/184

[56] References Cited

#### U.S. PATENT DOCUMENTS

Re. 34,507 1/1994 Egawa et al. . Re. 34,599 5/1994 Susyznski et al. . Re. 34,789 11/1994 Fraden .

(List continued on next page.)

#### FOREIGN PATENT DOCUMENTS

66021/86	of 1986	Austria .
778199	2/1972	Belgium
1258052	8/1989	Canada .
1265355	2/1990	Canada .
1314407	3/1993	Canada .

(List continued on next page.)

#### OTHER PUBLICATIONS

J. W. Moore et al., "Noncontact tympanic thermometer", *Medical & Biological Engineering & Computing*, vol. 16, No. 5, Sep. 1978, pp. 580–584.

D. E. Lees et al., "Noninvasive Determination of Core Temperature During Anesthesia", *Sourthern Medical Journal*, vol. 73, No. 10, Oct. 1980, pp. 1322–1324.

J. Fraden et al., "Application of Pyro–Electric Polymer Film to Medical Thermometry," *Proceedings of the Eighth Annual Conference of the IEEE/Engineering in Medicine and Biology Society*, 86CH2368.9, vol. 3 of 3, Forth Worth, Texas, Nov. 7–10, 1986.

J. Fraden, "Application of Piezo/Pyroelectric Films in Medical Transducers," *Journal of Chemical Engineering*, vol. 13, No. 3, Mar./Apr. 1988, pp. 133–138.


(List continued on next page.)

Primary Examiner—Christopher W. Fulton Attorney, Agent, or Firm—Nixon & Vanderhye P.C.

### [57] ABSTRACT

A disposable probe cover for an electronic tympanic ("ear") thermometer is made of closed-cell foam that is laminated with a sheet of high density polyethylene plastic film substantially transparent to infrared radiant energy. The probe cover is hollow, defining a central passageway that accepts an oversized probe and can potentially accommodate a variety of different probe configurations. The foam construction allows the probe cover to stretch over the oversized probe—retaining the probe cover on the probe during use. The foam construction also allows the probe cover to deform when inserted into a patient's ear-sealing the patient's ear canal to prevent external heat and light from affecting the measurement, and providing a high degree of cushioning and comfort not available in any prior disposable probe cover design. The film laminate and foam provide a barrier that is impervious to germs, fluids and other secretionsfurther reducing the already low risk of cross-contamination. The probe cover design facilitates accurate, repeatable measurements using a temperature measuring instrument based on receiving and analyzing the radiant infrared heat energy emitted by the tympanic membrane (eardrum).

### 8 Claims, 11 Drawing Sheets

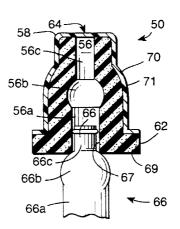

	U.S. PA	FENT DOCUMENTS	3,738,892	6/1973	Curcio .
D 210 051	0/1070	Cata	3,742,191		Poole et al
D. 218,851	9/1970	Everest.	3,750,471		Bremer .
D. 246,766 D. 254,959		Everest .	3,777,568		Risgin et al
D. 300,728	4/1989		3,781,748		Bishop et al
D. 303,008		O'Hara et al	3,781,837	12/1973	Anderson et al
D. 317,414		Lanci et al	3,798,366		Hunt et al
D. 318,812		Matsuura et al	3,809,228 3,809,229		Fowler et al Wahlig .
D. 321,487		Manno .	3,809,920		Cohen et al
D. 329,389		Hines .	3,812,847		Moore et al
D. 329,395	9/1992	Mackay .	3,822,593		Oudewaal .
D. 329,396	9/1992	Mackay .	3,822,598	7/1974	Brothers et al
D. 336,862		Ayton et al	3,832,669	8/1974	
D. 337,534	7/1993	Swift.	3,833,115	9/1974	1
D. 337,954		Makita et al	3,834,238	9/1974	
D. 338,412		Curbbun .	3,838,600		Ersek et al
D. 342,681		Mackay .	3,849,530	11/19/4	Wyeth et al Cornett, III et al
D. 370,860 738,960	6/1996 9/1903	Pompei et al Vaughan et al	3,851,029 3,878,836	4/1975	Twentier.
1,363,259	12/1920		3,880,282	4/1975	
2,696,117		Harrison .	3,929,018	12/1975	Turner .
2,804,069	8/1957	Schwamm et al	3,942,891	3/1976	Spielberger et al
2,844,031		Rosenthal .	3,949,740	4/1976	Twentier .
2,848,998		Bryan .	3,987,899	10/1976	Vyprachticky .
2,871,701	9/1959	Knudsen .	3,999,434	12/1976	Yen.
2,877,500	9/1959		3,999,537	12/1976	Noiles .
2,904,480		Rainer et al	4,005,605	2/1977	
2,969,141		Katzin .	4,022,855		Hamblen .
2,972,991		Burke .	4,024,397 4,054,057		Weiner .
3,023,398 3,054,397	2/1962	Siegert . Benzinger .	4,054,057	10/1977 12/1977	Kluge . Essen .
3,156,117	11/1964	Benzinger .	4,062,239	12/1977	Fowler et al
3,179,805	4/1965	Astheimer .	4,081,678	3/1978	Macall .
3,190,436	6/1965	Diamant .	4,091,922	5/1978	Egler .
3,193,978	7/1965	Bader .	4,117,926	10/1978	Turner et al
3,234,593	2/1966	Lerner et al	4,148,304	4/1979	Mull .
3,277,715	10/1966	Vanderschmidt .	4,159,766	7/1979	Kluge .
3,282,106		Barnes .	4,166,389	9/1979	Montren .
3,301,394		Baermann et al	4,166,454	9/1979	Meijer .
3,335,715		Hugenholtz et al	4,168,626	9/1979	Fullager .
3,349,896		Ensign et al Ensign et al	4,183,248 4,191,197	1/1980	
3,367,186 3,368,076	2/1968 2/1968	Clifford .	4,191,197	3/1980 3/1980	Benzinger . Wacker .
3,465,149	9/1969		4,197,944		Catlin .
3,469,449		Keller .	4,201,222		Haase .
3,469,685	9/1969		4,226,910	10/1980	
3,491,596	1/1970	Dean .	4,233,512	11/1980	Rupert .
3,500,280	3/1970	Ensign .	4,241,828	12/1980	Bourdelle et al
3,507,153		Jones et al	4,271,358		Schwarz.
3,526,135		Wortz .	4,275,591	6/1981	
3,531,642		Barnes et al	4,297,685		Brainard, II .
3,531,992		Moore . Wortz .	4,301,682		Everest .
3,581,570 3,605,750		Sheridan et al	4,312,357 4,315,150		Andersson et al Darringer et al
3,626,757		Benzinger .	4,341,992		Goldstein .
3,641,345		Coackley et al	4,343,182		Pompei .
3,650,153		Schwab .	4,343,185		Knute.
3,653,263		Poole et al	4,350,166		Mobarry .
3,663,917	5/1972	Mahmoodi .	4,351,616	9/1982	Farnstrom et al
3,673,868		Beury, III et al	4,362,166	12/1982	Furler et al
3,678,751		Mead et al	4,372,690		Berman et al
3,681,991		Eberly, Jr	4,378,489		Chabinsky et al.
3,701,347		Belkin .	4,379,971		Smith et al
3,703,892	11/1972	Meyers .	4,380,998	4/1983	Kieffer, III et al
3,719,396 3,724,448	3/1973 4/1973	VanDeWalker et al Lima .	4,392,005 4,400,341	7/1983 8/1983	Mohrman . Sorensen .
3,729,998		Mueller et al	4,400,341	11/1983	
3,735,864		Eckhart .	4,420,265		Everest et al
3,738,172	6/1973		4,425,921		Fujisaki et al
3,738,173	6/1973		4,433,924		Quinn, III .
3,738,479	6/1973	Sato .	4,436,438	3/1984	Voznick .

4,454,370	6/1984	Viznick .	5,088,834 2/1992 Howe et al
4,456,390		Junkert et al	5,094,544 3/1992 Ignatowicz.
4,457,633		Andrews .	5,127,742 7/1992 Fraden.
	9/1984		5,150,969 9/1992 Goldberg et al
4,471,354			5,153,563 10/1992 Goto et al
4,475,554		Hyndman .	5,159,936 11/1992 Yelderman et al
4,481,417	11/1984	Inglee .	
4,487,208	12/1984	Kamens .	5,163,418 11/1992 Fraden et al
4,493,564	1/1985	Epstein .	5,167,235 12/1992 Seacord et al
4,494,881	1/1985	Everest .	5,169,235 12/1992 Tominaga et al
4,509,522	4/1985	Manuccia et al	5,172,978 12/1992 Nomura et al
4,510,115	4/1985	Gokcen et al	5,178,464 1/1993 Fraden .
4,515,165	5/1985	Carroll .	5,179,936 1/1993 O'Hara et al
4,524,779		Brown, Jr	5,183,337 2/1993 Pompei.
4,527,896		Irani et al	5,188,459 2/1993 Mino et al
4,537,791		Tarjan .	5,199,436 4/1993 Pompei et al
4,566,808		Pompei et al	5,229,612 7/1993 Pompei et al
4,572,365		Bruno et al	5,232,284 8/1993 Egawa et al
4,588,306		Burger et al	5,264,375 11/1993 Bang et al
4,602,642		O'Hara .	5,271,407 12/1993 Pompei et al
			5,292,347 3/1994 Pompei.
4,607,963		Ulrickson .	5,293,862 3/1994 O'Hara et al
4,614,442		Poncy .	5,293,877 3/1994 O'Hara et al
4,619,271		Burger et al	5,313,951 5/1994 Zhao .
4,626,686		Pompei et al	5,319,202 6/1994 Pompei.
4,634,294		Christol et al	5,325,863 7/1994 Pompei .
4,636,091		Pompei et al	5,333,784 8/1994 Pompei .
4,644,163	2/1987	Selander.	
4,652,145	3/1987	Bjornberg .	
4,659,234	4/1987	Brouwer .	5,352,038 10/1994 Schmidt et al
4,662,360	5/1987	O'Hara et al	5,352,039 10/1994 Barral et al
4,679,949	7/1987	Sakamoto .	5,358,333 10/1994 Schmidt et al
4,684,018	8/1987	Jarund .	5,368,038 11/1994 Fraden .
4,691,712	9/1987	Brown, Jr	5,381,796 1/1995 Pompei.
4,727,500		Jackson et al	5,388,907 2/1995 Aoyama et al
4,763,522		Pompei .	5,391,001 2/1995 Rupert et al
4,765,752		Beynon et al	5,404,125 4/1995 Mori et al
4,784,149		Berman et al	5,411,032 5/1995 Esseff et al
4,790,324		O'Hara et al	5,445,158 8/1995 Pompei .
4,797,840		Fraden .	5,458,121 10/1995 Harada .
	1/1989		5,469,855 11/1995 Pompei et al
4,801,212			5,479,931 1/1996 Mooradian.
4,823,949	4/1989		5,487,607 1/1996 Makita et al
4,831,258		Paulk et al	5,515,847 5/1996 Braig et al
4,854,730		Fraden .	5,516,010 5/1996 O'Hara et al
4,859,079		Wickersheim et al	5,522,662 6/1996 Shiokawa .
4,863,281		Suszynski .	5,528,041 6/1996 Pompei .
4,874,253		Pompei et al	5,609,564 3/1997 Makita et al
4,895,164	1/1990	Wood.	
4,900,162	2/1990	Beckman et al	5,645,350 7/1997 Tang 128/664
4,907,895	3/1990	Everest .	FOREIGN PATENT DOCUMENTS
4,911,559	3/1990	Meyst et al	TOKEROW MILEWI BOCOMEWIS
4,914,673	4/1990	Imura .	0098402 1/1984 European Pat. Off
4,919,505		Bartosiak et al	0445783A2 9/1991 European Pat. Off
4,932,789		Egawa et al	0674162A2 9/1995 European Pat. Off
4,955,980		Masuo .	715359A1 6/1996 European Pat. Off
4,986,672		Beynon .	2167973 8/1973 France.
4,993,419		Pompei et al	2343234 9/1977 France.
4,993,424		Suszysnki et al	4422974A1 1/1995 Germany.
5,001,657		Yagura et al	•
		-	19604201A1 8/1997 Germany.
5,011,296		Bartosiak et al	55-154426 12/1980 Japan .
5,012,813		Pompei et al	56-161134 12/1981 Japan .
5,017,018		Iuchi et al	56-167428 12/1981 Japan .
5,017,019		Pompei .	57-35739 2/1982 Japan .
5,018,872		Suszynski et al	57-35740 2/1982 Japan .
5,024,533		Egawa et al	57-35741 2/1982 Japan .
5,031,619		Pompei .	57-212039 12/1982 Japan .
5,046,482	9/1991	Everest .	59-135439 9/1984 Japan .
5,051,590	9/1991	Kern et al	60-187829 9/1985 Japan .
5,051,595		Kern et al	100319 6/1987 Japan .
5,054,936	10/1991		147818 11/1962 U.S.S.R
5,056,682		Meyst et al	1425765 2/1976 United Kingdom .
5,066,142		DeFrank et al	1518521 7/1978 United Kingdom.
5,081,359		Pompei .	WO93/03666 3/1993 WIPO .
5,081,998		Yelderman et al	WO93/19662 10/1993 WIPO .
5,001,550	1/1/74	resormali et al	11 O25/12002 10/1223 WILO.

WO95/14913	6/1995	WIPO .
WO95/18961	7/1995	WIPO .
WO96/07877	3/1996	WIPO .
WO96/19938	7/1996	WIPO .

#### OTHER PUBLICATIONS

J.M. Looney, Jr. et al., "Ear Thermometry," *Medical Electronics*, Jun. 1989.

M. Benzinger et al., "Tympanic Clinical Temperature," presented at the Fifth Symposium on Temperature, Washington, D.C., Jun. 21–24, 1971, sponsored by the National Bureau of Standards, American Institute of Physics, and Instrument Society of America.

Information Sheet, Model 1M, Thermopile Detector, Dexter Research Center, 2 pages, Oct. 1980, Michigan.

Information Sheet, Model 2M, Thermopile Detector, Dexter Research Center, 2 pages, Oct. 1980, Michigan.

Information Sheet, Model M5, Thermopile Detector, Dexter Research Center, 2 pages, Oct. 1980, Michigan.

Information Sheet, Model DR26, Dual Element Thermopile Detector, Dexter Research Center, 3 pages, Oct. 1980, Michigan.

Information Sheet, Model DR46, Thermopile Detector, Dexter Research Center, 1 page, Oct. 1979, Michigan.

Information Sheet, Model 1010, Low Noise Amplifier, Dexter Research Center, 2 pages, Oct. 1980, Michigan.

Dexter Research Center, Inc., Product Price List—Domestic, 1–1–83.

Y. Houdas, et al., "Human Body Temperature, Its Measurement and Regulation," p. 83, Plemum Press, New York and London (Date unknown).

Advertisement, "Optical Calibration," Det Tronics, Sunnyvale, CA, *In Tech*, p. 48, Oct. 1987.

P. Gaudet, "Omega Tympanic Infrared Temperature Measuring Instrument, Product Specification," Exergen Corporation, pp. 1–3 (Date unknown).


ASTM Designation: E 1112–86 (Reapproved 1991), "Standard Specification for Electronic Thermometer for Intermittent Determination of Patient Temprature," pp 1–4 (Reprinted from the Annual Book of ASTM Standards, Philadelphia, PA).


"Electronic Thermometers," *Medical Electronics*, Jun. 1996, pp. 118–120.


Schieferdecker, J., et al., "Infrared thermopile sensors with high sensitivity and very low temperature coefficient," *Sensors and Actuators* A 46–47 (1995) 422–427 (printed in the Netherlands).

P.C. Lanchester, "Digital thermometer circuit for silicon diode sensors," *CRYOGENICS* 1989 vol. 29 Dec., (Received May 30, 1989), Southampton, UK, pp1156–1159.

Trial Exhibit list from Thermoscan, Inc. V. Sherwood Medical Co. (Date unknown).

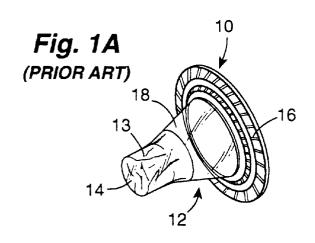

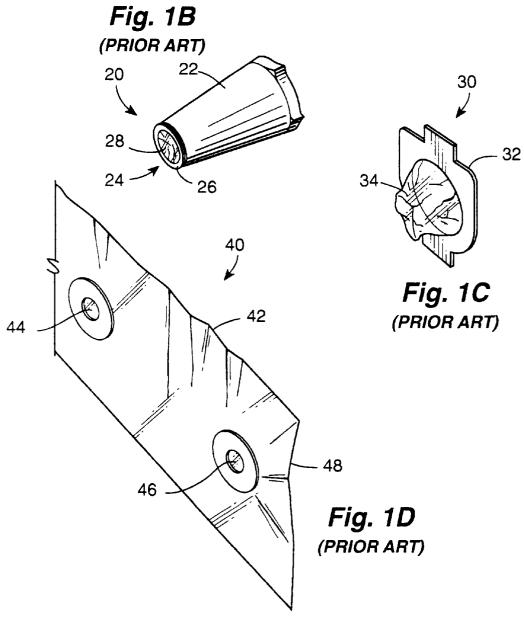

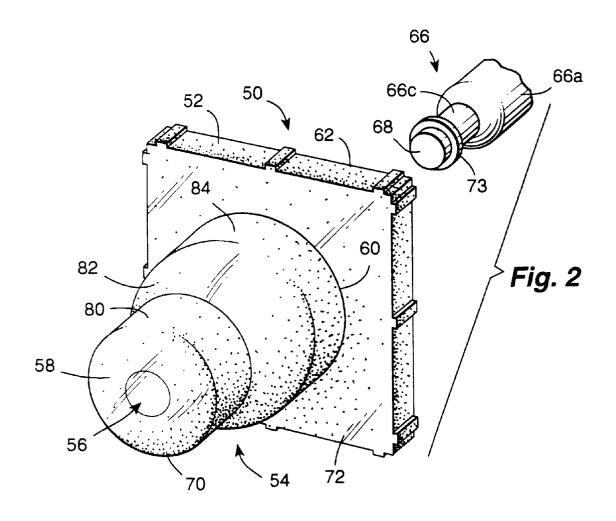

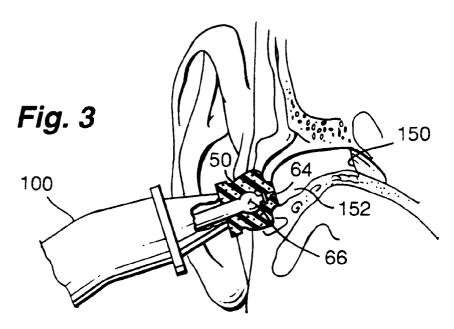
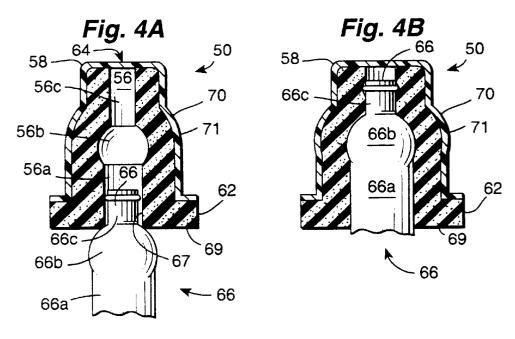


Fig. 3A 150




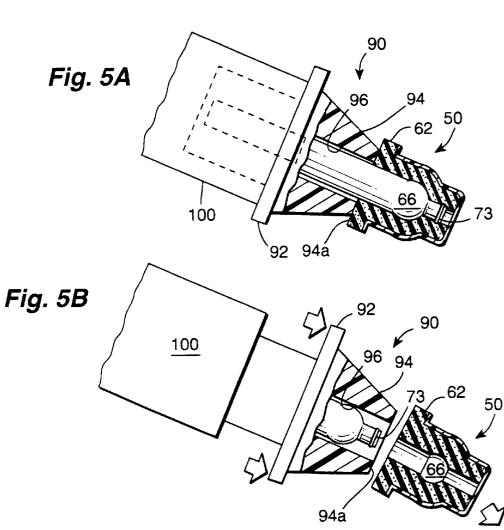



Fig. 6A

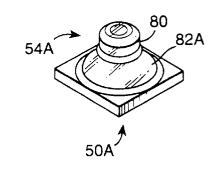

50F 54F 62F **Fig. 6F** 

**Fig. 6B**54B
80
86
82B

50B

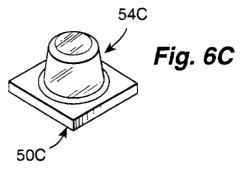

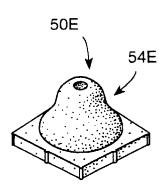
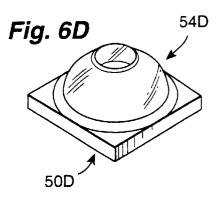
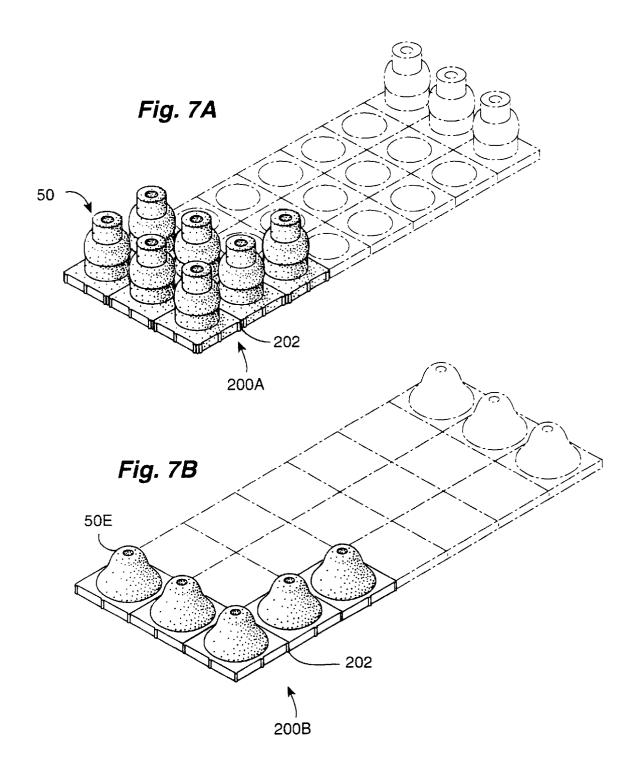


Fig. 6E









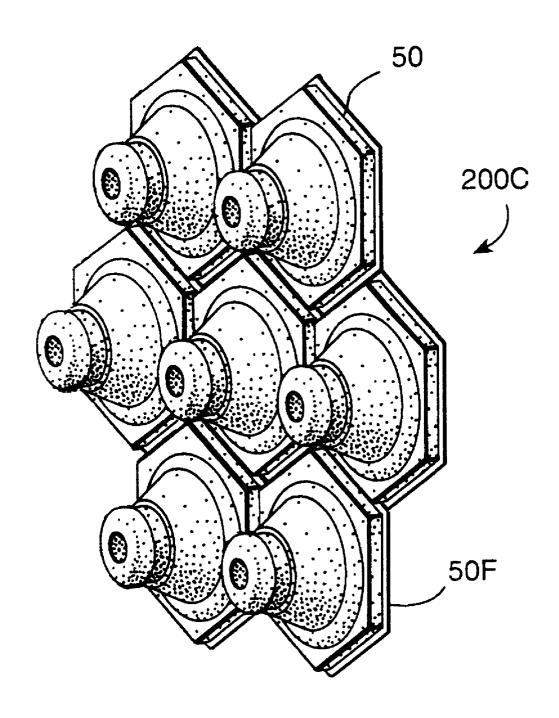


Fig. 7C

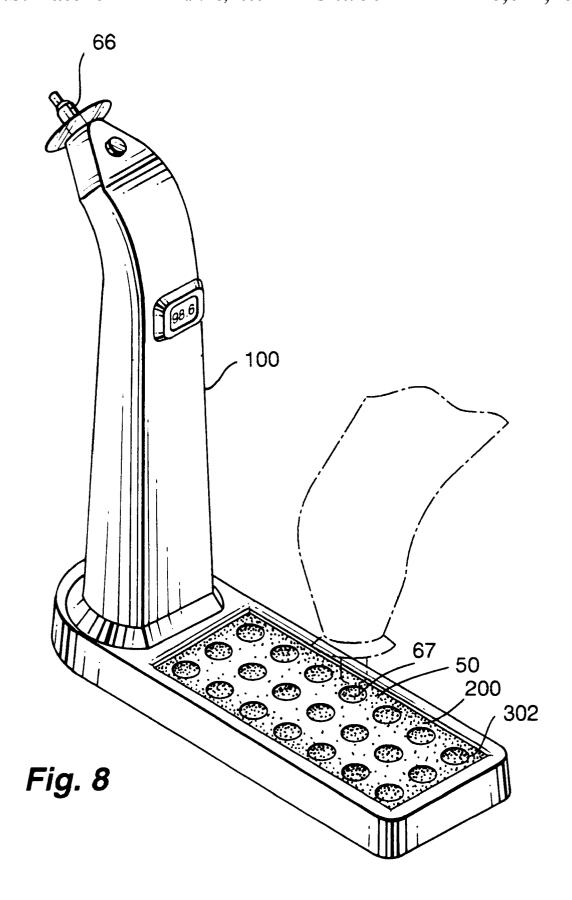

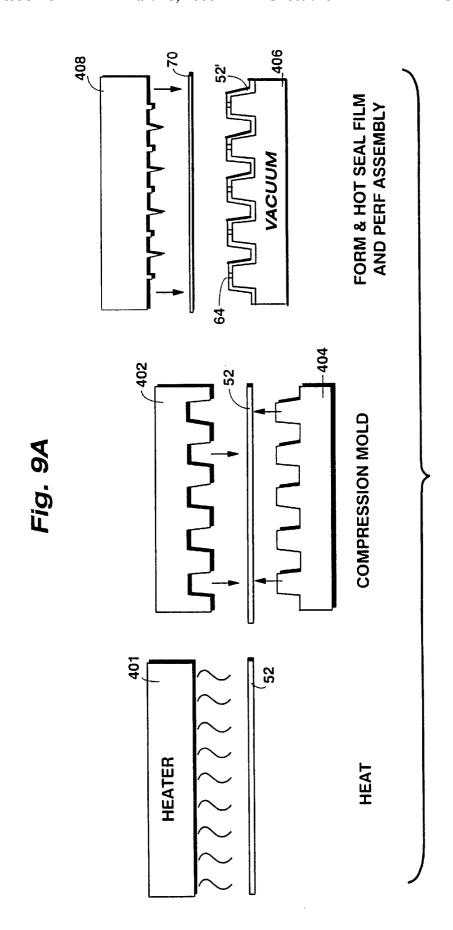

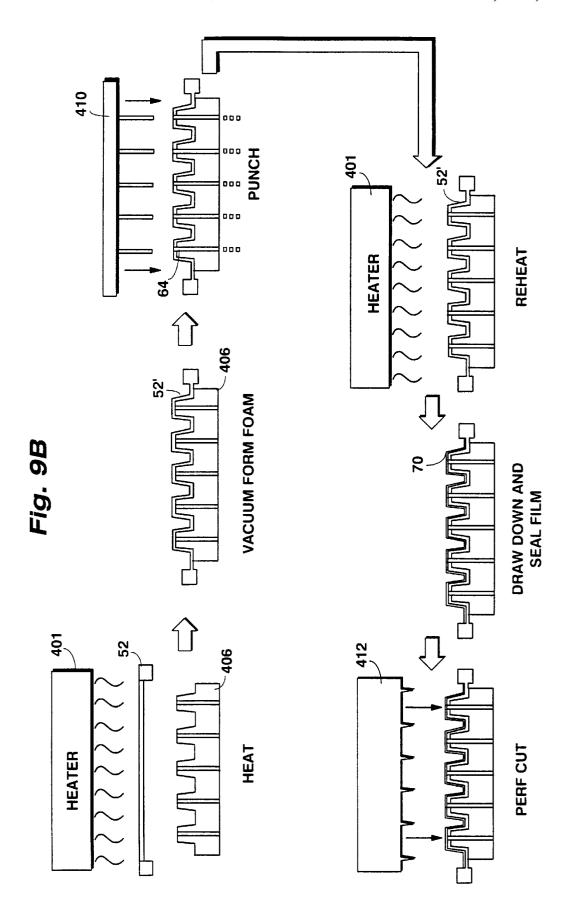

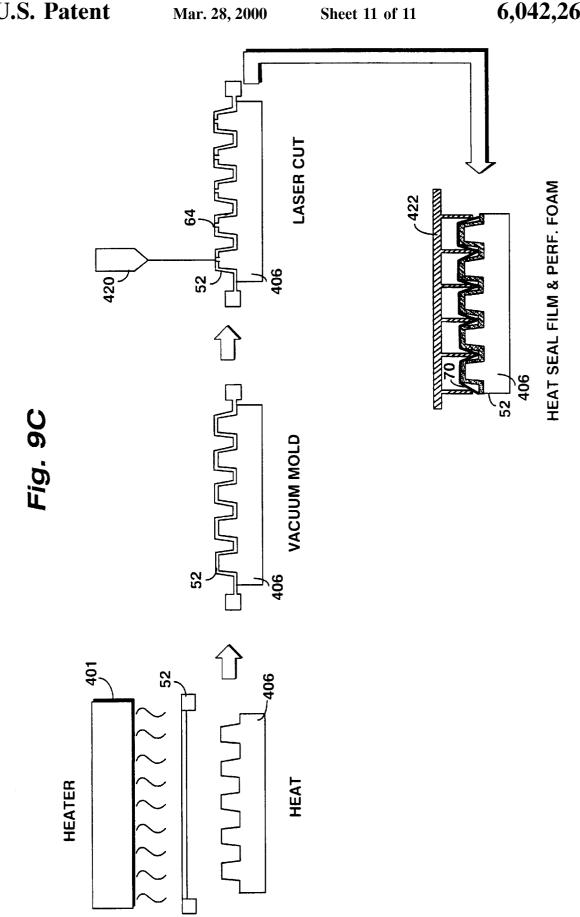

### TYMPANIC THERMOMETER PROBE **COVER**

This is a division of application Ser. No. 08/747,423, filed Nov. 12, 1996, now U.S. Pat. No. 5,833,367.

This invention relates to a probe cover for a thermometer that measures temperature by sensing the radiant energy given off by the eardrum ("tympanic membrane"). More particularly, the invention relates to a disposable foam probe of an ear thermometer temperature measuring system.

#### BACKGROUND AND SUMMARY OF THE INVENTION

Doctors, nurses, parents, and other care providers all need to be able to rapidly and accurately measure a person's body temperature. To find out whether a person is sick, the first thing a care provider usually does is take the person's temperature. Someone running a fever is likely to have an infection. A doctor or nurse can tell a lot about how a patient is doing by monitoring the patient's temperature over time and noting how it has changed.

There are three kinds of thermometers in wide use today: glass thermometers,

"electronic" thermometers, and

ear ("tympanic") thermometers.

Glass thermometers are very inexpensive, very small and easy to store, and don't require batteries or other special supplies. For this reason, glass thermometers are probably 30 the most widely used temperature measuring device in the home. However, glass thermometers have the disadvantage that they are very slow in making measurements—they typically require several minutes to reach body temperature. This is uncomfortable for the patient, and may be very 35 troublesome when it is necessary to take the temperature of a small child or an invalid. In addition, glass thermometers are typically accurate only to within a degree, may be susceptible to errors in placement, and can be broken easily.

Because of these disadvantages, most hospitals and doctors' offices now use instruments commonly known as "electronic" thermometers. Most of us have had our temperature taken by an electronic thermometer at one time or another. The electronic thermometer includes a portable, hand-held battery powered unit with a display, and a sepa- 45 pile" or pyroelectric heat sensor. This special heat sensor rate probe. A wire usually connects the probe to the handheld unit. The probe is long and thin, and has the same general shape as a glass thermometer. To use this kind of electronic thermometer, a nurse first covers the probe with a long thin disposable plastic probe cover that completely 50 covers the probe. The disposable probe cover helps prevent the spread of disease by avoiding direct contact between the reusable probe and the germs in the patient's mouth. The nurse then puts the end of the probe under the patient's tongue. An electronic temperature sensor within the probe 55 electrically senses the patient's temperature, and sends a signal to a microcomputer in the hand-held unit. The handheld unit usually beeps when the temperature measurement is finished, and displays the patient's temperature on the display. The nurse can then remove the probe from the patient's mouth, strip the probe cover off the probe, and throw away the used disposable probe cover.

This type of electronic thermometer has achieved wide acceptance in hospitals because it is reasonably accurate, can be used with familiar placement techniques, and is (because 65 of its disposable, replaceable probe covers) easily reusable for a number of different patients. Although the electronic

hand-held unit is itself more expensive than most households are willing to pay, the overall cost of using this kind of electronic thermometer is relatively low because the disposable probe covers are inexpensive (two to three cents per cover, for example) and a single hand-held electronic unit may last for years and can be used to take the temperatures of many thousands of patients.

Electronic thermometers offer speed, ease of reading, and accuracy improvements over glass thermometers, and also cover that increases the comfort, accuracy and repeatability 10 eliminate the possibility of mercury poisoning. Although such electronic thermometers have achieved a fair degree of success, they have certain significant disadvantages. For example, they need to be constantly calibrated, are relatively easily broken, and often require a relatively long time (thirty seconds or more in many cases) to make an accurate measurement. There are also problems with taking a temperature from the patient's mouth due to breathing, keeping the thermometer under the patient's tongue, etc. Crosscontamination of infectious diseases is also a concern 20 because the mouth is a "wet orifice."

> More recently, a new kind of electronic thermometer has appeared on the market. This new kind of thermometer works by measuring the temperature of your eardrum. Since the eardrum is also known as the "tympanic membrane," 25 these thermometers are sometimes called "tympanic thermometers."

Why the eardrum? The carotid artery that supplies blood to the hypothalamus—the body's temperature control center—passes through the eardrum. For this reason, the temperature of your eardrum corresponds very closely to the core temperature of your body. Although doctors and scientists have known this fact for many years, only since the mid-1980's have commercial devices been available to measure eardrum temperature in a clinical setting.

Ear or "tympanic" thermometers work by receiving and analyzing the radiant heat ("infrared") energy coming from the eardrum. Just as you can feel the heat when you hold your hands up in front of a warm fire, a tympanic thermometer can detect eardrum temperature without having to actually touch the eardrum by receiving the radiant heat energy coming from the eardrum.

Commercially available tympanic thermometers consist of a portable, hand-held battery powered main unit providing electronics, a display and a probe containing a "thermocalled a "thermopile" is especially sensitive to the eardrum's radiant heat energy. Microelectronics can determine eardrum temperature from the electrical signals provided by the "thermopile" sensor. The thermopile's sensing probe or "nacelle" is typically an integral part of the tympanic thermometer's main unit—reducing the potential for breakage of the sensor assembly and (at least potentially) increasing reliability and accuracy.

The sensing probe typically has a tapered or cone shape for easy insertion into the outer ear canal. All of the commercially available tympanic thermometers accept disposable probe covers. These disposable probe covers further minimize the risk of spreading disease, and also promote cleanliness by preventing ear wax and other secretions from 60 contacting the probe. Different manufacturers have different probe cover configurations (see FIGS. 1A-1B), but all of these probe covers have the common characteristic that they are relatively transparent to (and thus do not block or reduce) the radiant heat energy given off by the eardrum.

To use the ear thermometer, a nurse or other care provider inserts a disposable probe cover onto the instrument's sensing probe. Once the disposable probe cover is in place, the

nurse or other caregiver inserts the covered probe into the patient's outer ear and then presses a button to command the instrument to make a measurement. The measurement time is usually very rapid—on the order of two seconds or less. The patient's temperature instantly shows on the instrument's display. The instrument may then be removed from the patient's ear, and the disposable cover can be stripped off the instrument and discarded.

Ear thermometry has advantages over other temperature measuring techniques:

The measuring time is very rapid—usually less than two

Because the ear is a dry orifice, cross-contamination is not much of an issue-and individual, disposable probe covers further reduce the already low crosscontamination risks.

The theoretical accuracy of the measurement is very high (for example, on the order of one tenth of one degree). Because of the short measurement time and the use of the ear as the measuring point, it is possible to rapidly measure the temperature of children, invalids and sleeping patients—and in other situations where it is difficult to get a patient to sit still for thirty seconds with a probe under their tongue.

Despite these many clear advantages, ear thermometry has not yet achieved wide success in the medical market- 25 place. Even though many hospitals are believers in the concept of ear thermometry, the hospital market overall has converted less than twenty-five per cent of its temperature measurements to ear thermometry—and the hospitals that have converted are often displeased with their choice.

The main reason for past failures is that existing ear thermometer/probe cover combinations do not provide the high, repeatable accuracy required in a demanding hospital environment. Nurses are often unable to duplicate ear thertemperature twice with existing commercial ear thermometer/probe cover combinations, you may get two very different readings. Since accurate, repeatable temperature measurements are important or even critical to medical diagnosis and treatment (for example, to detect a 101.5° F. hospital fever threshold or to establish a temperature pattern over time), it is important for temperature measurements to be as accurate and repeatable as possible.

For an ear thermometer to make an accurate, repeatable measurement, the ear thermometer probe must be aimed 45 directly at the eardrum. Readings will be different unless the orientation of the probe within the ear is duplicated nearly exactly from one reading to the next. Because different structures within the ear give off, absorb and reflect different amounts and types of radiant energy, even slight differences 50 FIGS. 1A and 1B can help a nurse insert the probe into and in alignment between the probe sensor and the patient's eardrum can cause significant differences in the temperature

Poor technique is the leading cause of ear thermometer inaccuracy and non-repeatability—and disposable probe 55 cover design is an important factor in the nurse's ability to repeatably, correctly align the probe within the patient's ear. There have been several general approaches to tympanic thermometer probe cover design in the past. In one design approach, the disposable probe cover has a generally rigid structure designed to guide and place the probe (see FIGS. 1A & 1B for examples). In another design approach, the probe cover has very little structure, acting solely as a barrier against contamination and playing no role at all in probe placement (see FIGS. 1C & 1D for examples).

FIGS. 1A-1D show some examples of prior art tympanic thermometer probe covers. Prior art FIGS. 1A and 1B show

examples of prior art probe covers that are designed to assist in guiding and placing the probe. The FIG. 1A prior art probe cover 10 is made by Thermoscan, Inc. of San Diego Calif. and marked with "U.S. Pat. No. 5.088,834". This probe cover 10 is made out of a single unitary piece of thin, lightweight, semi-rigid, stiff hollow translucent plastic material shaped into a cone 12. The cone 12 includes a thicker "shank" portion 18, and terminates in a thinner end portion 13 comprising an integral piece of thin plastic film 14 that 10 is impervious to moisture but transparent to (and does not absorb) radiant heat energy at the infrared wavelengths emitted by the eardrum. The FIG. 1A one-piece probe cover 10 is designed to be inserted over a correspondingly coneshaped probe. A surrounding retaining ring base 16 retains 15 cover 10 on the probe during use. When cover 10 is inserted over the probe, plastic film 14 is stretched tightly over the end of the probe to provide a thin, wrinkle-free film "window" that is substantially transparent to infrared radiant energy. The conically shaped stiff "shank" portion 18 of cover 10 tapers to glide into the passageway of a person's ear leading to the eardrum—thus helping to position the probe relative to the eardrum. This stiff construction also allows probe covers 10 to be stacked one inside another for compact shipment.

FIG. 1B shows a different prior art probe cover 20 sold by Sherwood Medical Co. of St. Louis, Mo. (this design is also explained in Sherwood's U.S. Pat. Nos. 5,516,010, 5,293, 862, 5,179,936, 4,790,324, 4,662,360 and/or 4,602,642 and possibly also in a still-pending continuing application of these). The FIG. 1B probe cover 20 comprises a thick rigid plastic cone 22 for insertion into a person's ear. Probe cover 20 is made out of injection molded, rigid polyethylene or polypropylene white plastic. Small tabs or "ears" within the hollow inner space formed within plastic cone 22 retain mometer readings. If you try to measure the same person's 35 cover 20 on a specially designed corresponding probe. Plastic cone 22 terminates in a circular opening 24. A retaining ring 26 formed around opening 24 retains a thin transparent polypropylene or polyethylene plastic film membrane 28. Film membrane 28 is in an unstretched and wrinkled condition when the cover is not in position over a probe. When rigid probe cover 20 is inserted over a probe, the film membrane 28 stretches over the probe end to remove the wrinkles and provide a thin, substantially transparent film layer of uniform thickness. This stretched film layer 28 allows infrared radiant energy to pass through and reach the infrared probe inside rigid probe cover 20 without letting fluid, earwax, germs or anything else from the ear cavity come into contact with the probe.

> The rigidity of disposable probe covers 10, 20 shown in position the probe within the patient's ear. The rigidity of these probe covers 10, 20 also facilitates easy ejection of the probe covers from the probe after measurement. Unfortunately, this same rigidity makes these probe covers exceedingly uncomfortable for the patient. Because every person's ear has a different opening size, the "Lone size fits all" approach of the FIG. 1B probe cover 20 leads to discomfort and other problems. For example, probe cover 20 shown in FIG. 1B is known to become stuck in a person's ear, and can also scratch the inside of the ear. Frequently also, the probe tip will puncture the membrane film of the FIG. 1B probe cover when the clinician is inserting the thermometer into the probe cover.

The FIG. 1A probe cover 10 design is somewhat less 65 uncomfortable than the FIG. 1B design, but the rather large diameter of conical shank 18 can be uncomfortable for some people. Moreover, its unitary, one-piece construction

requires sophisticated techniques for drawing thin film portion 14 from the relatively thicker shank portion 18 during manufacture without puncturing the film. Also due to their rigidity, the FIG. 1A and FIG. 1B designs do not allow for total sealing of the ear canal to prevent external light and heat from entering during a measurement—degrading the accuracy and repeatability of temperature measurement.

FIGS. 1C and 1D show example prior art probe covers that follow a different approach of being mechanically transparent so as not to affect probe placement in any way. 10 The FIG. 1C probe cover 30 is manufactured by Diatek, Inc. of San Diego, Calif. See U.S. Pat. No. Re. 34,599. This probe cover 30 includes a flat paper backing 32 that supports a piece of transparent film 34. Film 34 is stretchable, and completely envelopes the probe when cover 30 is placed 15 over the probe-stretching tightly over the probe so as to provide a radiation-transparent film layer between the probe and the patient's ear. The FIG. 1C design has several advantages, including inexpensive construction, flat storage, and stackability. A big disadvantage of the FIG. 1C probe 20 cover 30 is that the probe cover is of no assistance in reliably, repeatably positioning the probe within the person's ear. This probe cover 30 is designed to be used with a temperature sensing "gun" also manufactured by Diatek (see U.S. Pat. No. 4,863,281)—an unpopular design because of 25 the bad connotations associated with putting a gun to someone's head.

FIG. 1D represents another prior art probe cover 40 manufactured by Exergen Corp. of Newton, Mass. and marked "U.S. Pat. No. 4,993,419." This disposable probe 30 cover 40 comprises a film sheet 42 that is dispensed in a pre-perforated roll. In use, the nurse inserts a pin portion of the ear thermometer positioned opposite the probe into a first hole 44, and then wraps film 42 around the front of the probe and continues wrapping until a second hold 46 is aligned 35 the nurse that the probe is properly in position. with the same pin portion. The nurse can then tear sheet 42 from the rest of the roll alone the line of perforations 48. Holes 44, 46 retain sheet 42 on the probe during measurement. Although probe cover 40 is inexpensive to manufacture, it is difficult for a nurse to position properly onto the probe. In addition, the FIG. 1D probe cover 40 may not always reliably wrap around the front portion of the probe—creating wrinkles that can effect measuring accuracy. The FIG. 1D probe cover 40 can also slip off the probe without the nurse realizing.

The FIG. 1D cover also has repeatability problems. In actual use, clinicians regard repeatability as implying accuracy. If the same reading can be obtained from the same ear, within several tenths of a degree, then the thermometer is regarded as accurate, regardless of the actual core tempera- 50 ture of the patient. All of the manufacturers currently recommend that if a second temperature measurement is required, that it should be made in the other ear. The reason for this is not obvious, but it does have a tremendous bearing on a major problem concerning accurate measurements.

Simply stated, because the nacelles of tympanic thermometers are generally at or near room temperature, which can be up to forty degrees below patient core temperature, the nacelle can rapidly draw off heat from the ear canal, impacting successive measurements dramatically. An excellent demonstration of this effect is to place the Exergen unit in a patient's ear and leave it in place for about a minute, while taking temperatures every ten seconds during this time frame. The average decrease in readings is always several degrees over a minute. The Exergen models produce the 65 a sheet of high density polyethylene. most dramatic thermal reductions because the heat transfer is taking place between the ear canal, a FIG. 1D thin layer

of film, and a large stainless steel mass. To a lesser extent, the FIG. 1A-1C probe cover designs exhibit the same, effect, albeit to a lesser extent.

As evidenced by the example prior art tympanic thermometer probe covers shown in FIGS. 1A-1D, there is a long felt but unsolved need for an ear thermometer disposable probe cover that:

is helpful to the nurse or other care provider in correctly, repeatably positioning the probe relative to the patient's eardrum.

is inexpensive to make,

is easy to place over the probe before use and easy to strip from the probe after use,

is comfortable for the patient, and

promotes measuring accuracy.

No prior solution meets all of these needs.

The present invention provides a new disposable probe cover for an ear thermometer that meets all of these needs and provides other advantages as well. The disposable probe cover provided by this invention is made from a compressible, deformable material such as plastic foam. The disposable probe cover provided by this invention has the following advantages:

Very easy to use.

Easy to strip off probe after use.

Much more comfortable to the patient.

Unlikely to lodge in the ear canal.

Foam provides a spring action to push itself out of the helping to prevent it from becoming lodged in the ear. Will not cause scratches or other abrasive injury.

Completely conformable, allowing a comfortable sealing of the ear canal.

The patient can detect proper sealing by the effect it has on the change of hearing—allowing the patient to confirm to

Prevents external light and heat sources from adding inaccuracies to the temperature measurement.

Provides excellent thermal isolation between the ear canal and the thermometer probe body-improving measurement accuracy by reducing thermal drift associated with the cold junction of a thermal pile-type measurement instrument.

Is an excellent thermal insulator and barrier. By limiting the amount of heat drawn off the ear canal, the ability to achieve repeatable measurements is greatly enhanced.

Provides extremely low emissivity in the infrared wavelength band of interest—further reducing error factors introduced through reflection within the ear canal.

Is inexpensive to manufacture.

Stretches over the probe to retain the probe cover without any special retention tabs or other mechanisms.

Can stretch to fit a variety of different configuration

Provides conforming resiliency that accommodate a 55 wide range of different ear sizes without being uncomfortable.

Can, in one embodiment, provide an externally-bonded thin-walled pre-stretched barrier (which can be less than 0.0075 inches) of high-density polyethylene that is substantially transparent to infrared radiation of interest while providing an industry accepted form of infection resistance and control.

Can be inexpensively vacuum-formed from a sheet of closed-cell foam that is dye cut and further laminated with

Allows use of a variety of different oversized probe configurations.

In one embodiment, preferred black or darker color foam provides greater repeatability by blocking extraneous light, absorbing stray ambient light, and having very low emissivity with negligible thermal mass.

#### BRIEF DESCRIPTION OF THE DRAWINGS

These and other features and advantages provided by the present invention will be better and more completely understood by referring to the following detailed description of preferred embodiments in conjunction with the drawings of which:

FIGS. 1A-1D show some examples of prior art disposable ear thermometer probe cover designs;

FIG. 2 shows an example disposable probe cover pro15 temperature of which is to be measured.

Vided in accordance with this invention;

FIG. 3 shows an example of how the FIG. 2 probe cover can be used with a tympanic thermometer to make an accurate, repeatable temperate measurement;

FIG. 3A shows a more detailed view of how the preferred <sup>20</sup> embodiment probe cover provided by this invention seals the ear passage while providing sufficient resilient conformal support to assist a nurse in repeatably positioning the infrared probe relative to the person's eardrum;

FIGS. 4A and 4B show how the preferred embodiment probe cover can accept and stretch over an oversized probe;

FIGS. 5A and 5B show how the preferred embodiment probe cover can be easily ejected from the probe;

FIGS. 6A-6F show examples of alternate preferred 30 embodiment probe cover configurations;

FIGS. 7A–7C show examples of how multiple probe covers can be packaged; and distributed in a common, preperforated block;

FIG. 8 shows an example ear thermometer and associated <sup>35</sup> probe cover dispenser tray designed to accept the probe cover block shown in FIGS. 7A–7C; and

FIGS. 9A–9C show example processes for manufacturing the preferred embodiment probe cover.

# DETAILED DESCRIPTION OF PRESENTLY PREFERRED EXAMPLE EMBODIMENTS

FIG. 2 shows an example disposable probe cover 50 provided by this invention. Probe cover 50 provides a cover or sheath for covering an ear thermometer radiant energy sensing probe 66. Probe cover 50 prevents germs and secretions within a person's ear canal 152 (see FIG. 3) from contacting probe 66 while being substantially invisible to infrared radiant heat energy emitted by the person's eardrum that is directed along the path toward probe sensing end 68.

FIG. 2 shows probe 66 before it has been inserted into and covered by probe cover 50. FIGS. 3, 3A, 4A, 4B and 5A show probe cover 50 in cross-section after the probe 66 has been inserted into and covered by the probe cover.

Probe cover **50** in the preferred embodiment is made of closed-cell black plastic foam (e.g., a closed cell polyolefin foam such as polyethylene or polypropylene, thermoset closed-cell polyurethane foam, or expanded polystyrene foam). The foam may have a density of anywhere from 2 pounds to 6 pounds per cubic foot for example. In this embodiment, the foam wall thickness of the finished probe cover **50** will vary, but in the preferred embodiment may have a thickness of between ½6 and ¾6 of an inch to provide a desired cushioning and sealing effect. The foam construction of probe cover **50** allows the probe cover to deform when inserted into a patient's ear (see FIGS. **3** and **3**A)—

8

sealing the patient's ear canal and providing a high degree of cushioning and comfort not available in any prior disposable probe cover design. The closed-cell foam construction probe cover 50 also provides a barrier that is impervious to germs, fluids and other secretions—further reducing the already low risk of cross-contamination.

Probe cover **50** in this embodiment is hollow, defining a central passageway **56** that accepts oversized sensor probes **66** of potentially different configurations. The foam construction of probe cover **50** allows the probe cover to stretch over the oversized probe **66**—retaining the probe cover on the probe during use. The foam construction of probe cover **50** also provides a heat insulating effect—insulating the probe **66** from heat sources other than the eardrum the temperature of which is to be measured

In more detail, probe cover 50 a base portion 62 and a hollow tubular projection 54—which may both be molded from the same foam sheet 52. Base portion 62 may be square in shape (and, in one embodiment, may have dimensions of 0.80 inches by 0.80 inches). Hollow tubular projection 54 projects outwardly from square base portion 62 (so that, in one embodiment, the total length of probe cover 50 may be 0.88 inches). Hollow tubular projection 54 has two ends 58, 60. The first end 58 (closest to the viewer in FIG. 2) is for insertion into a person's outer ear canal 152—see FIG. 3. The other end 60 joins and is integral with base 62.

Tubular projection 54 includes a central hollow passageway 56 disposed axially within and extending its length. This passageway 56 extends the entire length of probe cover 50, and in one example is defined along the central axis of projection 54. The passageway 56 terminates at an opening 64 at the projection first end 58, and at an opening 67 defined in an outer surface 69 of square base portion 62 (see FIG. 4A).

In use, a radiant energy temperature sensing probe 66 is inserted into passageway 56 through the base portion opening 67. The projection first end 58 is then inserted into the person's ear canal 152 (see FIG. 3). Infrared radiant heat energy emitted by the person's eardrum 150 (see FIG. 3A) passes through opening 64 and impinges upon the probe's sensing end 68 (which is now encased within the probe cover 50).

In this example, passageway 56 is dimensioned and shaped to accept an oversized heat sensing probe 66. Central passageway 56 terminates at end 58 in a substantially circular opening or window 64. The sensor end 68 of probe 66 (which may also be circular—and, in one example, forms the end of a circular infrared waveguide) is aligned with and aimed at probe cover circular opening 64 when probe 66 is inserted into probe cover 50. Radiant heat (infrared) energy can pass through opening 64 to impinge on probe sensing end 68.

In this example, the probe cover outer surface 71 has a cushioning obstruction 82 that prevents the probe cover 50 (and thus the probe 66) from being inserted too far into ear canal 152 and also providing a relatively fixed spacing between the probe end 68 and eardrum 150 each time the probe is inserted into the ear canal 152. In more detail, the preferred embodiment probe cover tubular projection 54 may comprise three sections: an end, insertion section 80; a middle, obstructing section 82; and a base joining section 84. In this particular example, sections 80, 82, 84 may have different diameters and shapes.

For example, insertion section **80** may be substantially circularly cylindrical with a uniform diameter alone its length extending 0.13 inches radially outwardly from the

edge of circular window 64. Insertion section 80 is preferably shaped and dimensioned to conveniently and easily slide into ear canal 152 without discomfort.

Middle obstructing 82 in this example preferably has a diameter that is substantially larger than that of end section 80 so it won't slide into the ear canal 152. Middle, obstructing section 82 acts as a "stop" to limit how much of probe cover 50 that can be inserted into the ear canal 152. In this example, middle section 82 may have a spherical shape providing an abrupt increasing diameter such that it acts as 10 a stop—abutting against the ear lobe and/or other part of the auricle (see FIG. 3A) to limit penetration of probe cover 50 into ear canal 152—and also providing a sealing effect (discussed below). The rapidly increasing diameter of the probe cover middle, obstructing section 82 helps ensure 15 repeatable temperature measurements by positioning probe end 68 the same distance from eardrum 150 each time the covered probe 66 is inserted into a given patient's ear. This constant distance between probe end 68 and ear drum 150 is time.

Base joining section 84 in this example may have a diameter that is slightly reduced relative to the maximum diameter of middle section 82. The purpose of base joining section 84 is simply to join the rest of projection 54 to base

Although the particular shape and outer contours of probe cover 50 shown in FIG. 2 provides various advantages, modifications to this shape may be desirable in certain 30 situations. FIGS. 6A–6E show some possible alternate configurations for probe cover 50. These alternate configurations are preferably each made from a foam sheet but may provide different outer surface shapes. For example:

eliminates tubular base portion 84, and provides a cylindrical end portion 80 and a middle portion 82A shaped as a hemisphere.

FIG. 6B shows another example foam probe cover 50b end section 80 and tapered portion 82A. This flat ledge 86 may abut against the ear lobe or other part of the auricle to provide a positive stop preventing further insertion of probe cover 50 into ear canal 152.

The FIG. 6C probe cover 50c design provides an <sup>45</sup> extending projection 54c that is uniformly linearly tapered.

The probe cover 50d design shown in FIG. 6D provides a projection 54d that is tapered as an ellipsoid to provide a hemispherical shape.

FIG. 6E shows a still additional probe cover design **50***e* in which the projection 54e has a complex tapered shape.

FIG. 6F shows a probe cover design 50f in which projection 54f has a conical shape terminating in a bulbous end portion, and has a hexagon shaped base portion 62f.

In one embodiment, probe cover internal central passageway 56 may be formed to match the outer surface contour of probe 66. This feature is shown most clearly in FIGS. 4A & 4B. In the particular example shown for illustrative purposes, probe 66 comprises a cylindrical portion 66a, a bulbous center portion 66b, and a terminal cylindrical portion 66c that terminates in probe sensing end 68. In this particular example, the diameter of terminal cylindrical portion 66c is less than the diameter of cylindrical portion **66***a*—and the bulbous center portion **66***b* has a generally hemispherical shape having the same diameter as that of cylindrical portion 66a. The shape of passageway 56 gen10

erally matches the shape of probe 66 except that the outer dimensions of the probe are greater than the inside dimensions of the passageway. Thus, in this example the probe cover internal passageway 56 may similarly define a lower cylindrical hollow 56a, a bulbous center hollow 56b, and an upper cylindrical hollow 56c communicating with one another to form a continuous passageway. In this embodiment, the probe 66 is "oversized" with respect to passageway 56 in this example in that:

- (a) the inner diameter of probe cover internal passageway lower cylindrical hollow 56a is less than the diameter of probe cylindrical portion 66a;
- (b) the inner diameter or probe cover passageway bulbous hollow 66b is less than the outer diameter of probe bulbous portion 66b; and
- (c) the inner diameter of probe cover internal passageway upper cylindrical hollow 56c is less than the outer diameter of probe cylindrical terminal portion 66c.

When probe 66 is fully inserted into the probe cover helpful in reliably establishing a temperature pattern over 20 internal passageway 56, the probe cover 50 stretches and conforms around probe 66 (due to its foam construction) to form a tight friction fit. Compare FIG. 4A with FIG. 4B. In this particular embodiment, the probe bulbous portion 66b mates with the probe cover bulbous hollow 56b to hold the probe cover 50 on the probe 66. More specifically, in this particular embodiment, the foam material defining the probe cover bulbous hollow 56b expands and stretches around the probe bulbous portion 66b to prevent axial movement of the probe 66 relative to the probe cover 50 (and thus ear canal 152 and eardrum 150) during measurement. The tight friction fit provided by the stretching of the foam material over probe 66 is aided by a high surface tension between the foam and the probe 66 to provide good retention. The coefficient of friction between the foam surface and the probe 66 is high FIG. 6A shows a modified probe cover 50 a that  $_{35}$  due to the molecular interaction between the two surfaces. The stretchable characteristic of foam also cover **50** allows the cover the adapt to a variety of differently dimensioned (and/or shaped) probes 66.

In one particular example, the probe cover end 58 may be that provides a flat ledge 86 in the transition region between 40 covered by a sanitary barrier 70. Barrier 70 may, in one embodiment, cover the entire outer surface 72 of probe cover 50—or it can be placed across opening 64. Barrier 70 prevents germs, fluids, ear wax and other contaminants from passing through opening 64 into contact with probe 66. Barrier 70 in one embodiment may comprise a thin clear high-density pre-stretched polyethylene plastic film sheet. For example, barrier 70 may be a plastic film membrane in the range of 0.0005" and 0.010" thick. Barrier 70 is substantially transparent to infrared radiant energy in the wavelength range of 8 to 14 microns emitted by the human eardrum—allowing such radiant energy to pass substantially unimpeded from the eardrum through circular opening 64 and passageway 56 to strike probe sensing end 68.

In use, barrier 70 across circular opening 64 should be of uniform thin thickness without any bumps, wrinkles, distortions or other obstructions. In this particular example, barrier 70 comprises a plastic film that is pre-stretched across opening 64 in an entirely unwrinkled state at time of manufacture. Barrier 70 therefore has no wrinkles in it at any time. Since probe cylindrical portion 66c may in one example be oversized (and/or include an oversized portion 73) with respect to probe cover internal passageway cylindrical hollow 56c as discussed above, inserting the probe 66 into cover 50 may cause the barrier 70 to further stretch across opening 64—from an already pre-stretched condition to a further stretched condition. The expansion of probe cover 50 around oversized probe 66 (e.g., structure 73) may

cause the normally-stretched barrier 70 to stretch even further across window 64 to provide a uniformly thin (e.g., less than 0.0075") high-density barrier that is impervious to germs and other contaminants but substantially transparent to at least a part of infrared radiant energy within the wavelength range emitted by eardrum 150. Deformation of foam based probe cover 50 upon contact with the ear may cause barrier 70 to stretch still further. More specifically, when the probe cover 50 is inserted into ear canal 152, the cover deforms axially outwardly to further stretch the barrier 10 70. This additional stretching of barrier 70 does not have any adverse effects so long as the barrier doesn't stretch so much that it punctures. Since the thickness of barrier 70 during measurement is a factor in the temperature reading, it is desirable to have different individual probe covers 50 pro- 15 vide about the same amount of stretching for insertion of a given probe 66.

Probe sensing end 68 may be positioned any desired distance from barrier 70 when the probe 66 is fully inserted into probe cover **50**. However, since the barrier **70** can never 20 be a 100% transparent to radiant energy (and will generally exhibit some reflectivity characteristics), it is desirable to control the distance between the probe sensing end 68 and the barrier 70 so as to provide a higher degree of measurement repeatability. Therefore, in this particular example, the mating of probe bulbous portion 66b with the corresponding stretch-fit probe cover passageway bulbous hollow 56b causes probe sensing end 68 to be disposed some fixed axial, repeatable distance from (without contacting the barrier 70 pre-stretched across and covering the opening 64.

FIGS. 3 and 3A show how probe cover 50 can be used with a tympanic thermometry instrument 100 to make a temperature measurement. Probe 66 is first fully inserted into the probe cover central passageway 56 by inserting probe sensing end 68 into base opening 67 and then exerting 35 pressure. Because of its deformable, stretchable foam construction, probe cover 50 stretches over probe 66 to provide a friction fit that securely retains the probe cover on the probe's outer surface—and also compensates for differences in shape and dimensions for differently configured 40 probe heads.

Conforming the central passageway 56's inner surface to the probe 66's outer surface also exactly aligns the probe sensing end 68 with the probe cover window 64.

50 insertion end 58 into the external acoustic meatus (ear canal) 152 of the patient. Probe cover insertion end 58 easily slides into ear canal 152. Because of its compressible foam construction, probe cover 50 conforms to the shape and dimensions of the ear canal 152—provide a high degree of 50 patient comfort and also sealing the ear canal to prevent external light, air currents and radiant heat from entering the ear canal during temperature measurement. Such external heat, air currents and/or light can affect the temperature measurement. Additionally, since the foam has an extremely 55 low emissivity in the wavelength bands of interest, it further reduces error factors introduced through reflection within ear canal 152. For example, black foam provides very low emissivity and thermal mass, and also absorbs ambient light. White colored foam can also be used. The foam color can be chosen to provide optimal radiant and optical characteristics (i.e., passivity, emissivity, absorption and reflectivity).

As mentioned above, the foam construction also provides a useful cushioning effect that keeps the hard surfaces of probe 66 away from the ear canal 152, provides a generally 65 comfortable feel, and easily accommodates ear canals 152 having different sizes and shapes. See FIG. 3A. It keeps

patient discomfort to a minimum, and is unlikely to lodge in the ear canal 152 or cause abrasive injury. Anyone who has ever had their temperature taken with a prior art rigid probe cover 50 as shown in FIGS. 1A and 1B will remember how uncomfortable the experience was. The soft, supple foam construction of preferred embodiment probe cover **50** provides an altogether different, much more pleasing experience. With probe cover 50 inserted the patient may feel pressure and will also have his or her hearing muffled due to the probe cover's sealing effect (e.g., as if you put your ear up to a sea shell). Thus, the patient can tell when the probe cover 50 has been fully inserted—an additional benefit if the patient is taking his or her own temperature (or to help train nurses in using proper temperature measuring technique). The inserted probe cover **50** does not feel like a too-big, hard object—the sensation most people experience with the prior art probe cover designs of FIGS. 1A and 1B.

The preferred embodiment probe cover 50 provided by the present invention also provides a useful positioning function guiding and centering probe 66 within ear canal 152 to aim probe 66 more or less directly at eardrum 150. This centering and aiming function provided by probe cover 50 improves the measuring accuracy and repeatability of tympanic thermometer 100 by ensuring that probe 66 is aimed at eardrum 150 and assumes more or less the same position relative to the eardrum and the ear canal 152 each time the probe is inserted into the same patient's ear canal. While some of the aiming and guiding responsibility still falls upon the nurse's technique (for example, how he or she grasps the instrument 100 relative to the patient's ear canal 152), probe cover 50 is helpful in this guiding and placement function.

Probe cover 50 also performs an important insulating function. Because foam is a natural insulator, it provides excellent thermal isolation of the ear canal 152 and the probe 66 (in particular, the so-called "cold junction mass" within the thermopile sensor of the thermometer 100 nacelle). The largest degree of error in any thermopile-type measurement instrument is usually the thermal drift associated with the cold junction. The foam construction of probe cover 50 dramatically reduces the  $\Delta \text{Temp}/\Delta \text{Time}$  relative to contact with the ear canal 152-reducing such thermal drift and improving measurement accuracy and repeatability. The ability to achieve repeatable measurements from the same or The nurse or other care provider then inserts probe cover 45 different ear is greatly enhanced by the insulative effect of probe cover 50 limiting the amount of heat drawn off the ear canal.

> Once the measurement is completed, disposable cover **50** can be stripped from probe 66 and discarded. One way to strip cover 50 from probe 66 is simply to grasp the cover with the hands and pull it off the probe. However, it is desirable to make the stripping operation as simple, fast and easy as possible while minimizing contact between the nurse's hands and any portion of cover 50 that has been in contact with the patient. FIGS. 5A and 5B show how probe cover 50 can be easily stripped from instrument 100 after use through use of an ejecting ring structure 90.

> In this example, instrument 100 may be equipped with a retractable ejection ring structure 90. Ejection ring structure 90 may be moved between two positions: a retracted position (as shown in FIG. 5A) and an extended, ejecting position (as shown in FIG. 5B). Ejection ring structure 90 may include an annular grasping portion 92 and a conical stripping ring 94. Probe cover base 62 may abut an upper surface 94a of stripping ring 94 and the probe 66 may pass through a central cylindrical passage 96 within the ejection ring structure 90 when the probe cover 50 is fully inserted

onto probe 66 and the ejecting ring structure is in the retracted position. To strip probe cover 50 from probe 66, the user simply moves ejection ring structure 90 from its retracted position to its ejecting position. In the ejecting position, ejection ring structure 90 may extend fully over probe 66 (see FIG. 5B)—overcoming the friction between the probe cover and the probe and stripping off the probe cover.

It also is desirable to help the nurse insert probe 66 into a probe cover **50** as easily and rapidly as possible without 10 requiring the nurse to touch, with his or her hands, any part of the cover that will be placed into direct contact with the patient. FIG. 7A-7C show how multiple foam probe covers 50 can be distributed together in the form of a block 200 for use with a dispensing tray 300 shown in FIG. 8. By way of example, the block 200a shown in FIG. 7A comprises twenty-one individual probe covers 50 of the type shown in FIG. 2. FIG. 7B shows a block 200b of twenty-one probe covers 50e of the type shown in FIG. 6E. FIG. 7C shows a block **200***c* of probe covers **50**f of the type shown in FIG. **6**F, the hexagonal base portions 54f providing a very compact packing density so many probe covers can be provided in a small flat surface area. In each case, the individual probe covers 50 are attached together by small foam fingers extending across a cut or perforated line 202. Probe covers 50 can be easily separated from one another by application of a small amount of mechanical pulling force-but stay attached to one another during shipment and handling of block 200 as a whole.

FIG. 8 shows an example probe cover dispenser tray 300 30 designed to accept and dispense a probe cover block 200. In this example, the probe cover block 200 is inserted into a rectangular tray opening 302 defined by tray 300, with the probe cover projecting portions 54 facing downward into the tray. This leaves opening 67 into which probe 66 is to be 35 inserted completely exposed. To attach an individual probe cover 50 to probe 66, a user may grasp the tympanic thermometer instrument 100 and insert probe 66 into probe cover opening 67 (see phantom lines in FIG. 8). The nurse then pushes probe 66 downwardly into probe cover passage 56—deforming and stretching the probe cover as shown in FIG. 4B to provide a tight frictional, high surface tension fit between the probe cover and the probe.

The probe 66 may stretch the foam composition of probe cover 50—causing the already stretched barrier 70 to 45 procedures. They take your temperature many times a day in become further stretched. The probe head expanded portion 73 helps to further stretch barrier 70 without requiring any contact between the probe end 68 and the barrier 70. Gently pulling the thermometer 100 away from tray 300 will provide enough force to separate the individual probe cover 50 50 from block 200—leaving probe 66 fully inserted into probe cover passage 54. The advantage of using the dispensing tray 300 shown in FIG. 8 is that it allows a nurse to very rapidly cover probe 66 with a probe cover 50 without having to ever touch or otherwise contact the outer probe 55 with what is presently considered to be the most practical cover surface 71 that is to come into contact with the patient. As discussed above, the barrier 70 may become still further stretched from its already-stretched position upon insertion of probe cover 50 into the ear canal 152 due to the outward deformation of the probe cover upon contact with the 60 and scope of the appended claims. structures of the ear.

FIGS. 9A-9C show different example processes for manufacturing probe covers 50. In the FIG. 9A example, a foam sheet 52 of desired thickness (for example, between ½4<sup>th</sup> of an inch and ¼ of an inch—and preferably between 65 1/16 and 3/16 of an inch thick) is heated by a heater 401, and is then placed between two mating compression halves 402,

404. The compression mold halves 402, 404 are shaped to provide the desired probe cover shape and inner passageway 54 shape (see FIGS. 6A-6E). In this example, the compression mold process may be used to pinch form the passageway 54. The compression mold halves 402, 404 may the be opened and the resulting formed sheet 52' may then be held in place on a workpiece 406 (or on lower mold half 404) by a vacuum. A film sheet 70 of desired thickness is then laminated on top of the formed foam sheet 52' using a hot die 408 that also punches opening 64 and provides perforations 202 between individual probe covers 50—or a numerically controlled ("NC") laser cutting device can be used to cut the perforations between the individual probe covers. Film sheet 70 may be laminated over the entire molded foam sheet 52', 15 or it may also be cut so as to cover a smaller surface area—but at least all of openings 64.

14

In the FIG. 9B example manufacturing process, a foam sheet 52 is heated by a heater 401 and is then placed onto a heated vacuum workpiece 406 that vacuum forms the foam sheet into the desired shape. A punch 410 cuts openings 64 into the sheet 52 to form individual probe covers 50. The same or different heater 401 is used to reheat the formed, punched foam sheet 52'. A film sheet 70 is then drawn down and sealed to provide a laminated film seal covering openings 64. A perforator 412 is then used to cut the perforations between the individual finished probe covers 50 to allow for separation of the individual pieces.

In the FIG. 9C example, heater 401 is used to heat the foam sheet 52 and a vacuum workpiece 406 forms the foam sheet into the desired shape as shown in FIG. 9B. Alaser 420 is then used to cut the holes 64 (and can also be used to cut the perforation between probe covers). Film 70 is then placed over the foam 52 and a heated die 422 is lowered into position. Pressure exerted by heated die 422 pre-stretches film 70 over holes 64, and the heat from the die bonds film 70 to foam 52 to maintain the film in a pre-stretched, wrinkle-free state. The film 70 is not pre-heated except locally along the edges of the perforations. Because the film 70 is heated only locally, it does not sag or become loose or wrinkled in the neighborhood of hole 64 but instead remains tightly stretched across the hole.

When you are sick, you want to make sure your doctor has the most accurate information concerning your conditionand you also dread having to put up with uncomfortable the hospital. The last thing you want is to feel uncomfortable every time they take your temperature. The prior art ear thermometer probe covers (see for example FIGS. 1A and 1B) are extremely uncomfortable because of their rigidity and other shortcomings. They can even get stuck in your ear, and they can scratch your ear. The present invention provides a new foam probe cover design that is much more comfortable and effective as compared to prior approaches.

While the invention has been described in connection and preferred embodiment, it is to be understood that the invention is not to be limited to the disclosed embodiment, but on the contrary, is intended to cover various modifications and equivalent arrangements included within the spirit

We claim:

- 1. A disposable probe cover comprising:
- a compressible, deformable material formed to provide a compressible, deformable hollow body that is dimensioned to be at least partially inserted into the outer ear canal of a human being, said deformable hollow body having a passage defined therethrough, the passage

being dimensioned for, in use, accepting at least a portion of a radiant energy sensing structure so that said deformable hollow body stretches and thereby frictionally engages said portion of said radiant energy sensing structure while being readily removable therefrom after 5 use.

- 2. A disposable probe cover as in claim 1 wherein said compressible, deformable material comprises foam.
- 3. A disposable probe cover as in claim 1 wherein said probe cover when inserted into the outer ear canal, deforms 10 probe cover, when inserted into the external acoustical so as to substantially seal the outer ear canal.
- 4. A disposable probe cover as in claim 1 wherein the passage terminates in an opening, said opening admitting electromagnetic energy into the passage to strike said radiant energy sensing structure.
- 5. A disposable probe cover comprising a hollow body comprising a compressible, deformable material, the probe

cover, in use, stretching to accept and thereby frictionally engaging an oversized radiant energy sensing probe and for insertion into the external acoustical meatus of the human body, said probe cover being readily strippable from said energy sensing probe after use.

16

- 6. A disposable probe cover as in claim 5 wherein said compressible, deformable material comprises foam.
- 7. A disposable probe cover as in claim 5 wherein said meatus, deforms so as to substantially seal the external acoustical meatus.
- 8. A disposable probe cover as in claim 7 further including a thin plastic film membrane laminated onto at least a 15 portion of the body and pre-stretched across the opening.

## Adverse Decision In Interference

No. 6,042,266, Edward P. Cheslock, Eric L. Canfield, Richard K. Harris, TYMPANIC THERMOMETER PROBE COVER, Interference No. 104,708, final judgment adverse to the patentees rendered November 26, 2002, as to claims 1-7.

(Official Gazette January 14, 2003)